

19 -80

Ultrasonic Scaler Ultrasonic Periodontal Therapy System



More comfortable periodontal treatment

Heating system

Upgraded water path heating system, built-in water heater, automatic preheating. Warmwater make patients comfortable.

One click to clean

Sandblasting air path cleaning and scaling water path cleaning, one click to activate. Efficiently clear the residual powder or liquid, avoid clogging and corrosion.


7-inch Large Touch Screen

Upgraded large touch screen, more intuitive, easier to operate.

Voice broadcasting system

Add voice broadcast function, more intelligent.



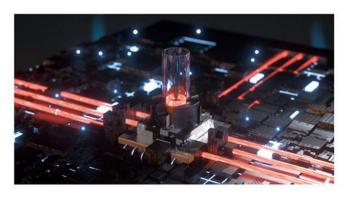

Large Touch Screen






One click to clean



Voice broadcasting system



Intelligent Multi-Core Chips

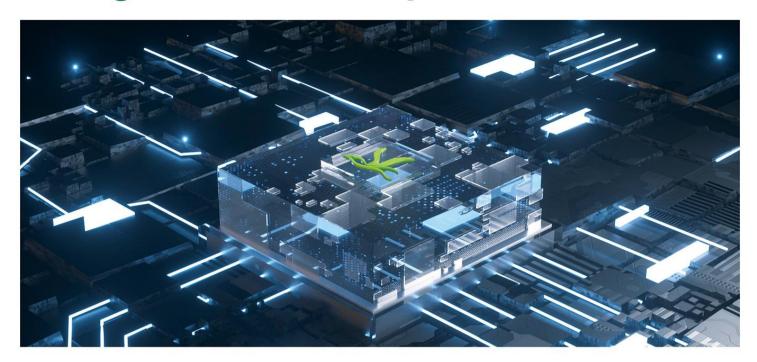

Built-in 32-bit ARM core processor, quick response, accurate feedback and intelligent adjustment of scaling power, automatic memory of latest setting.








Implant maintenance

It can effectively prevent and treat peri-implantitis, and protect the implant from damage

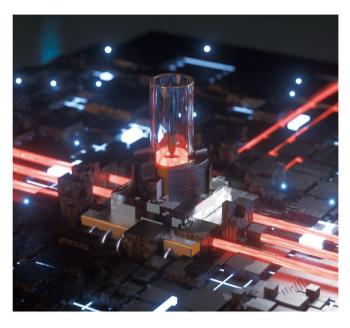

Heating system, more comfortable

Built-in water heater, automatic preheating.



Multifunctional wireless foot pedal

Wireless control, infinite possibilities, unconstraint

Effictive dental plaque control




Supragingival Air-Polishing Handpiece

0.7mm aperture air-polishing nozzle is made of steel (imported from Germany). Reliable and durable, air polishing is more concentrated, more efficient, more comfortable.

Regulate the powder flow under constant pressure

In the case of fixed pressure value, the doctor can scientifi-

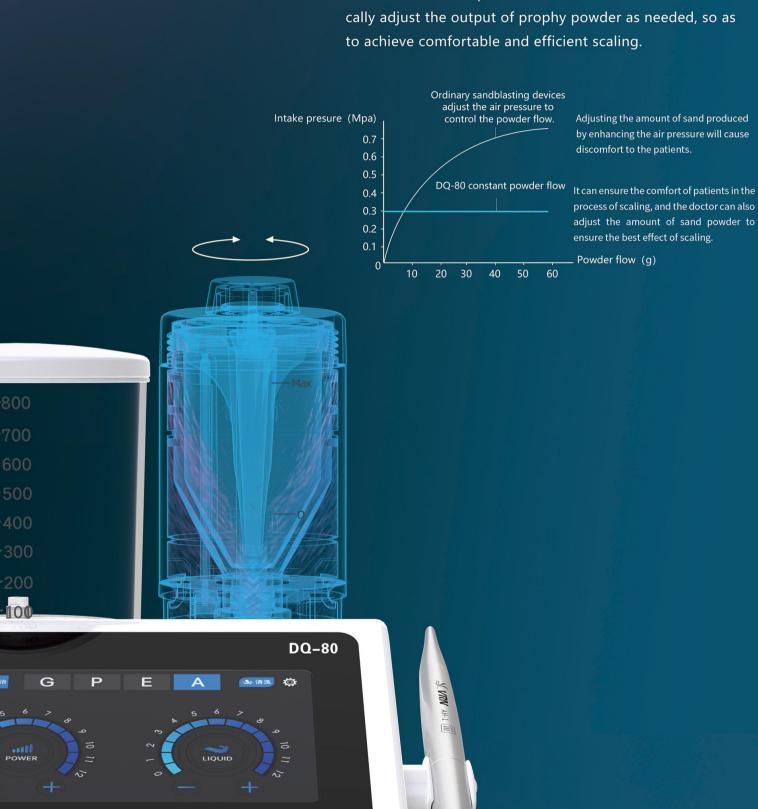

Comfortable ultrasonic scaling

The elliptical vibration track can realize the operation of 0° angle between the scaler tip and the tooth surface, which has the polishing effect.

It can effectively reduce the collision between the working tip and the teeth, achieve comfortable tooth cleaning.



Titanium Alloy Scaler Tips

Some scaler tips are made of titanium alloy material, which have lower hardness than enamel, no damage to the teeth surface, more comfortable.



Air polishing good partner



Supragingival Sodium Bicarbonate Supargingival Glycine Sub + Supragingival Glycine

SPECIFICATION:

DQ-80 Ultrasonic Scaler

Input Voltage: 230 V ~ or 120V ~, 50/60 Hz 0.5A

Main unit output: 25V~, 50Hz 2.8A

Input Power: 100VA

Tip amplitude: Minimum:1 μm, deviation-50%.

Maximum:100 µm, deviation+50%

Output half-excursion force: Minimum: 0.1 N, deviation-50%.

Maximum: 5 N, deviation + 50%

Tip vibrating frequency: 30 kHz ± 5 KHz

Tip output power: 3 W∼20 W

Air inlet pressure: 5 bar~6 bar (0.5 MPa~0.6 MPa)

Water temperature: 0°C~45°C

Main unit size: 310mm×370mm×200mm

Item	Qty.	Mark
Main Unit	1	DQ-80
Power Adaptor	1	
Wireless Foot Swtich	1	
Wired Foot Switch	1	
AH-1 Sandblasting Handpiece	2	
AH-2 Sandblasting Handpiece	1	
Subgingival Sandblasting Nozzle	40	
PH-2 Handpiece	2	
Water Tank	1	0.5
Needle	1	
Needle	1	
Waterproof Sticker	1	
Disinfection Box	2	
Sand Tank Components	2	
Self-locking Quick Connector	1	
Endo Wrench	1	
Wrench	2	0.5m
PU Pipe 4*6	2	2.5m
Braided Hose	1	
Self-locking Quick Female Joint	2	
Self-locking Quick Male Joint	1	
T type Quick Connector	2	PEG8-6
Fuse	1	
O Ring	10	11.1x1.8
	5 2	5x1 50x3
	6	8x1.5
	6	6.5x1.5
	6 8	11x1.5
Tips	20	P11x2 P12x2 P12Lx2 P12Rx P16x2 G1x1 G5x1 P1x1 IM1x1 E1x1
Instruction Manual	1	
QC Certificate	1	
Tips Instruction Manual	1	
Scaler Tips abrasion Contrast Card	1	
The Specification of O Ring	i	
Sand Dredging Components	1	
Endo Files	6	25#

Guilin Veirun Medical Technology Co., Ltd.

Guilin Veirun Medical Technology Co., Ltd. Add.: D-07 Information Industrial Park, National High-Tech Zone, Guilin Guangxi 541004 P. R. China

Tel.: +86-773-2260519 Email: info@veirun.com www.veirun.com



@VRNDENTAL Scan to visit VRN public page